


Digital Scanning Analyzers for Screening of Sexually Transmitted Diseases (STD)



Vision A1

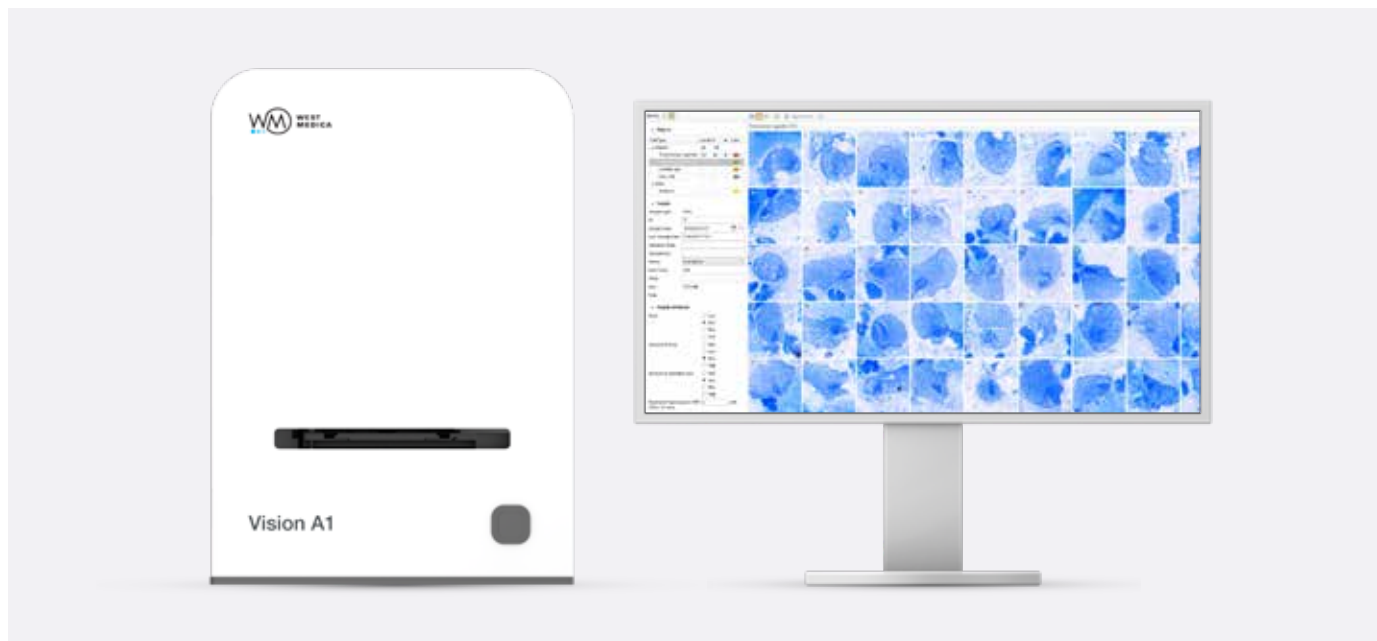

Key features

- Loading 1 slide
- Automatic sample scanning
- Automatic analysis and pre-classification
- For small and medium-sized laboratories

Scanning process




01

Place the slide to the slide holder



02

Install into the slide scanner and start scanning

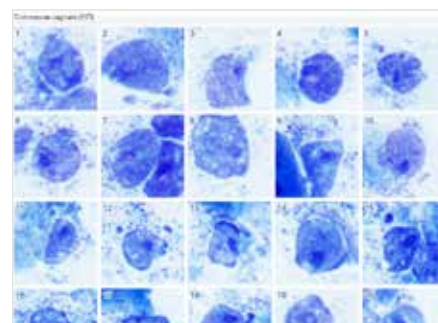

03

Get the result!

Vision Assist



- Loading 1 or 4 slides
- Automatic sample scanning
- Automatic analysis and pre-classification
- For small and medium-sized laboratories

Vision Pro




- Loading 4 or 8 slides
- Automatic sample scanning
- Automatic analysis and pre-classification
- Random access
- For medium and large laboratories

Vision Ultimate



- Loading 120 or 200 slides
- Automatic sample scanning
- Automatic analysis and pre-classification
- High scanning speed and image capture of the whole sample
- Random access and testing of urgent samples
- Night shift work without an operator
- For centralized and large-sized laboratories

Vision Cyto STD (Sexually Transmitted Diseases) Automatic screening for STD diagnostics

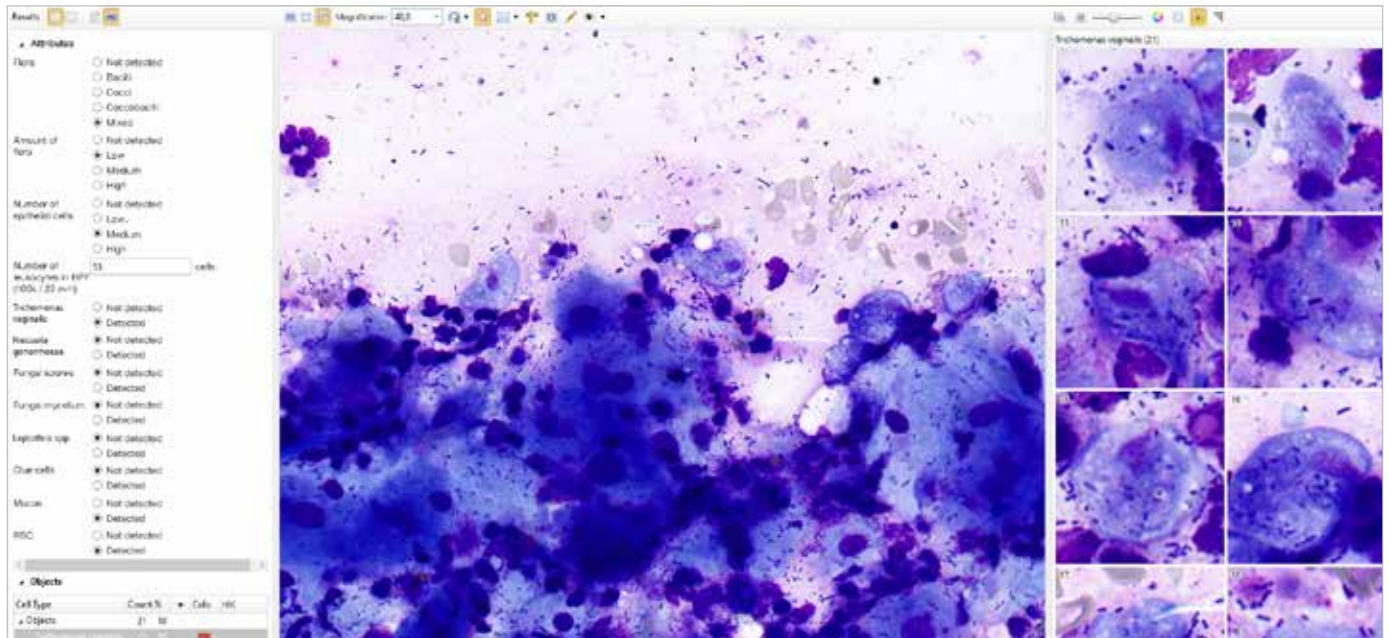

Key features

- Automatic sample scanning
- Creation of digital slides
- Testing for sexually transmitted infections

Automatic pre-classification


- Trichomonas vaginalis
- Neisseria gonorrhoeae
- Fungal spores
- Fungal mycelium
- Leptothrix spp.

Features of Vision Cyto STD



“Hotkeys”

Verification is done with one click

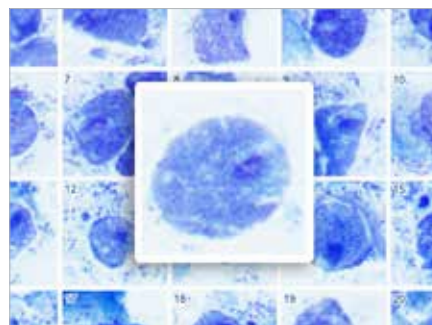


Quick view

Review any particular cell with higher magnification, fast evaluation of morphology


Fungal mycelium	Not detected		
Leptothrix spp.	Not detected		
Clue cells	Not detected		
Mucus	Not detected		
RBC	Not detected		
+ Objects			
Cell Type	Count	%	Reference Range, %
Objects	238	100	
Trichomonas vaginalis	238	100	
Neisseria gonorrhoeae	--	--	
Fungal spores	--	--	
Fungal mycelium	--	--	
Leptothrix spp.	--	--	
Clue cells	--	--	
Other	--	--	
Artifact	--	--	

Test results


Provided in a form of a table with information about identified cells and composition of microbiota

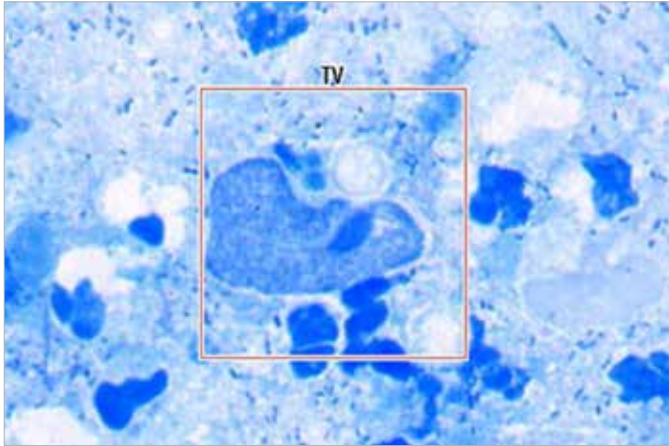


Trichomonas vaginalis

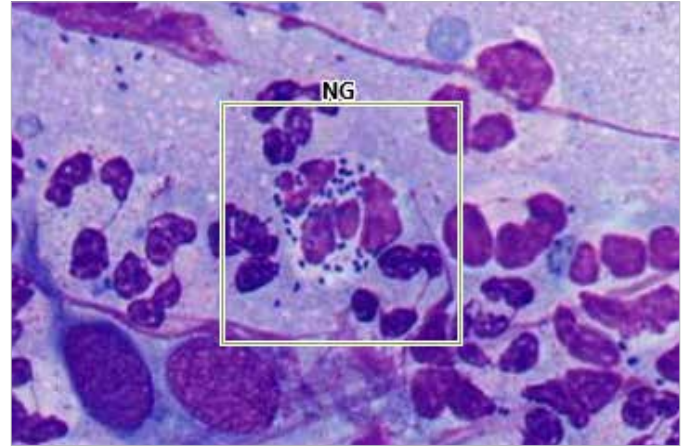


Neisseria gonorrhoeae

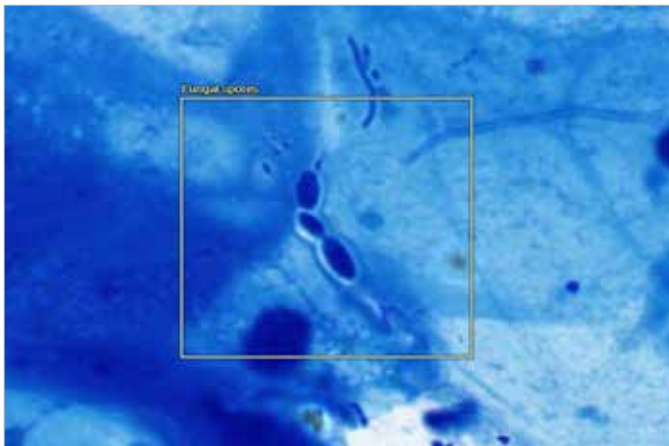

Fungal spores

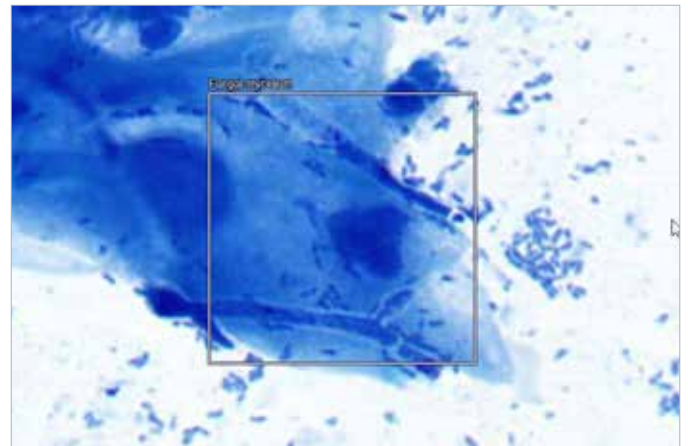


Fungal mycelium

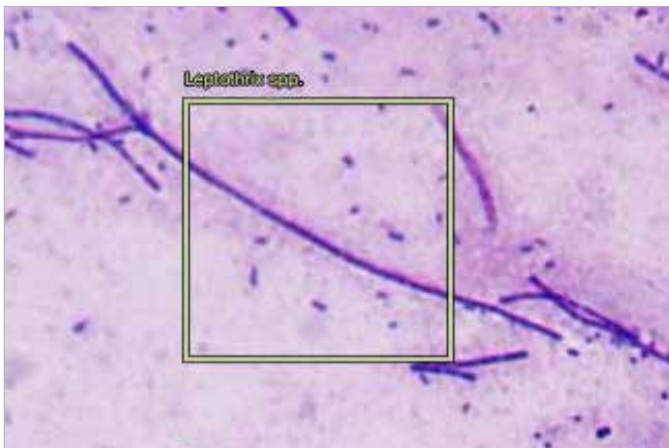

Leptothrix spp.


Vision Manager

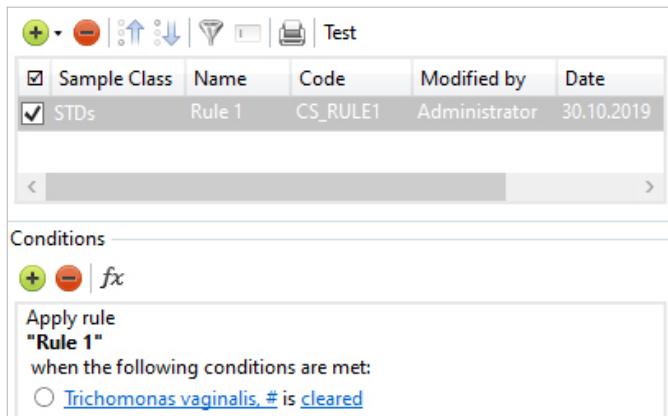

www.wm-vision.com

Data processing rules for analysis procedure automation


Vision Remote

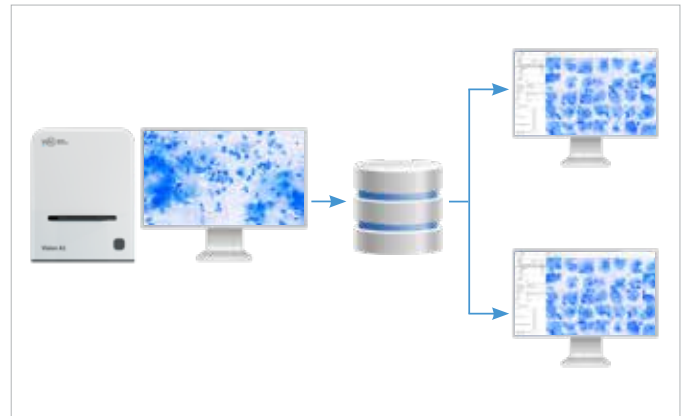

www.wm-vision.com

Remote workstation

Vision Suite

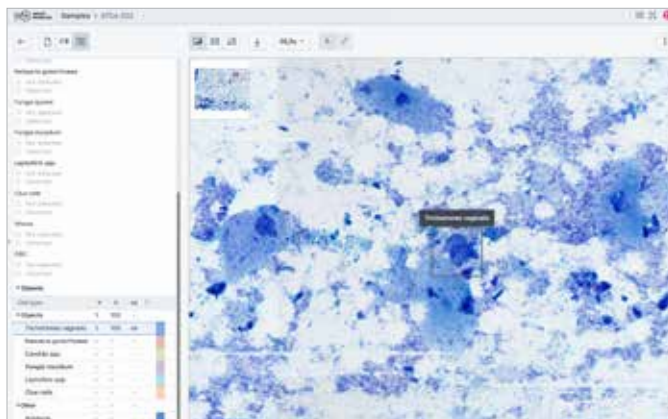

www.vision-suite.com

Remote operation of Vision software via Internet browser

Vision Expertise



www.vision-expertise.com


Online testing and validation of specialists' knowledge in microscopy

Artificial intelligence for digital microscopy



With artificial intelligence technologies, Vision's digital microscopy solutions optimize and standardize laboratory workflow.

Pathoview

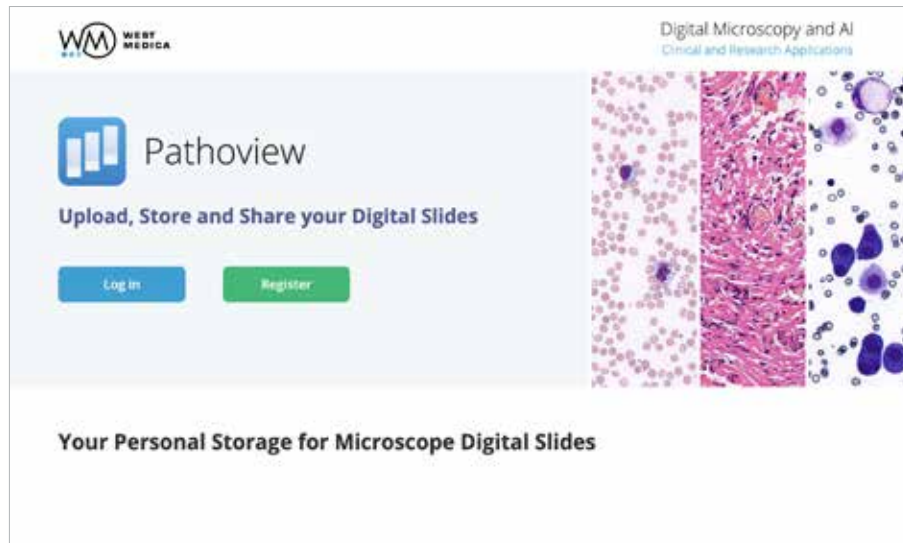

www.pathoview.com

Working with digital slides via internet browser

Upload, view, store and share digital slides



Uploading

Upload your slides in any convenient way: via a web browser, from the Vision and Vision Slide Viewer applications, or send them using your Vision Suite organization.



Storage

Save digital slides in Pathoview, and they will always be at your fingertips wherever you are. Ensure access to digital slides from any device.



Shared access


Information exchange — share and discuss clinical cases with colleagues. Quick and easy sharing of links to slides of any size.




Digital slides

Support various digital slide formats: .mrxs, .svs, .tiff, .ndpi, .dcm.


Vision and Pathoview integration

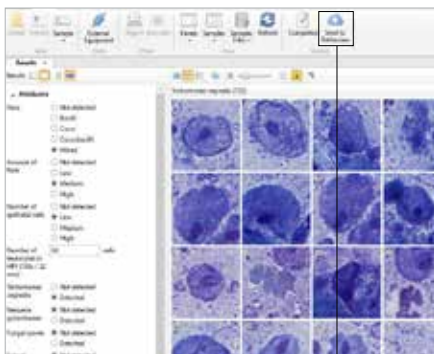

In Vision select a sample and click the “Submit” button on the toolbar

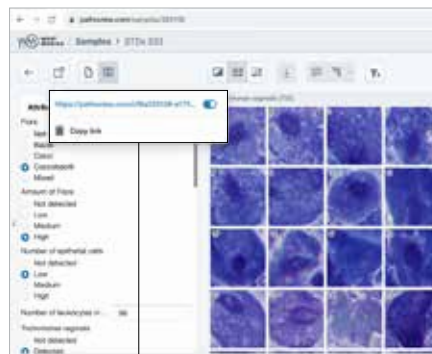

In Pathoview click the “Submit” button, copy the link and send it to your colleagues for consultation

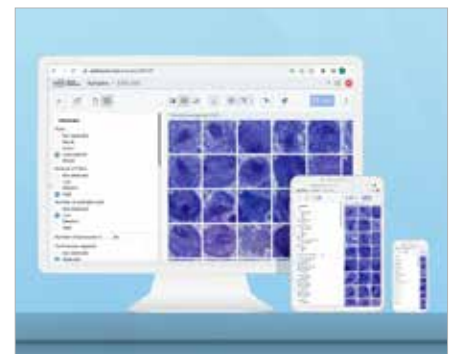

Your colleague views the digital slide on any device

Digital scanning analyzers

Key features	Vision A1	Vision Assist	Vision Pro	Vision Ultimate
Application modules	Vision Cyto STD	Vision Cyto STD	Vision Cyto STD	Vision Cyto STD
Working mode	Continuous loading	Queue (for 4 slides only)	Queue, random access	Queue, continuous and random access, STAT sampling, 24/7
Scanning	Automatic	Automatic	Automatic	Automatic
Loading a slide	1 slide	1 or 4 slides	4 or 8 slides	120 or 200 slides
Specifications	Scanning unit Personal computer Monitor	Scanning microscope Personal computer Monitor	Scanning microscope Personal computer Monitor	Scanning microscope Personal computer Monitor Control touch-monitor
Objectives	60x Oil	10x, 50x Oil, 100x Oil	10x, 40x/60x Oil, 100x Oil	2,5x, 10x, 63x/100x Oil
Research method	Bright field	Bright field	Bright field	Bright field
Illumination	Köhler illumination, LED	Köhler illumination, LED	Köhler illumination, LED	Köhler illumination, LED
Communications	Bidirectional connection to LIS, LIS2-A2 (ASTM), HL7	Bidirectional connection to LIS, LIS2-A2 (ASTM), HL7	Bidirectional connection to LIS, LIS2-A2 (ASTM), HL7	Bidirectional connection to LIS, LIS2-A2 (ASTM), HL7

Pre-classification has to be validated by the specialist reviewing the galleries of cells.

Please check the intended use and the registration in your local regulation.

Depending on the regulations of the respective region, some products may be used for research purpose only.

We reserve the right to change specification without notice.



West Medica Produktions- und Handels-GmbH
Brown-Boveri-Straße 6, B17-1
2351 Wiener Neudorf, Austria
tel.: +43 (0) 2236 892465, fax: +43 (0) 2236 892464
vienna@westmedica.com, www.wm-vision.com

Official Distributor